

ZEISS Optical Inserts

Sign up on MyZEISS Vision to

- Activate your free one-year warranty extension.*
- Track your personal eye health information and expand your vision knowledge.
- Save your unique pairing code online, as featured on the printed card.

*Valid if you register your ZEISS Optical Inserts within 90 days of purchase.



INSTRUCTIONS FOR USE
ZEISS Optical Inserts *Prescription*

The following instructions contain information on how to install and use ZEISS Optical Inserts *Prescription* hereinafter "ZEISS Optical Inserts" with Apple Vision Pro.

Read the instructions carefully before use.
Keep this document for your reference.

Apple Vision Pro User Guide:
<https://support.apple.com/guide/apple-vision-pro>

Contents

1. Legend

2. Package contents

3. Intended use

4. Safety notes

5. Installation and use

6. Storage

7. Maintenance and cleaning
8. Troubleshooting

9. Technical information


10. Warranty

11. Packaging label







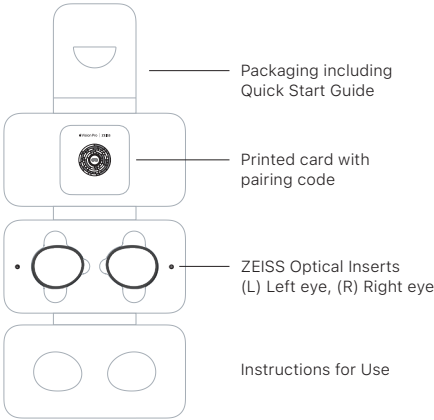






12. Disposal and recycling

13. Contact information

2 Package contents



1 Legend

 WARNING Death or serious injury Type and source of danger and possible harm. • Solution to avoid danger, etc.		
 WARNING Injury Type and source of danger and possible harm. • Solution to avoid danger, etc.		
 CAUTION Potential risk or product damage Type and source of danger and harm or possible damage. • Solution to avoid damage, etc.		
 Manufacturer's trade name and address (includes legal manufacturer/ specification developer)		
 Manufacturing date	 Medical Device	 Catalog number
 CE Mark (Certification)	 Unique device identifier	 Serial Number
 Identification of packaging material	 Magnetic field	 Pacemaker

3 Intended use

ZEISS Optical Inserts are intended to provide vision correction for a better user experience while wearing Apple Vision Pro. Users who normally wear eyeglass lenses or who otherwise require vision correction, can use ZEISS Optical Inserts as directed in this document to improve their visual experience.

ZEISS Optical Inserts are designed for use solely with Apple Vision Pro. ZEISS Optical Inserts cannot be used while driving a vehicle, operating heavy machinery, or in any other situation that requires attention to safety.

ZEISS Optical Inserts are not a substitute for prescription eyeglass lenses, reading eyeglass lenses, contact lenses, or any other vision correction eyewear that is intended for normal wear and use outside of Apple Vision Pro.

Never use the ZEISS Optical Inserts for vision correction outside of Apple Vision Pro, as such use could lead to unsafe situations. View the Apple Vision Pro Instructions for Use at: <https://support.apple.com/guide/apple-vision-pro>

4 Safety notes

 **WARNING | Death or serious injury**

1. ZEISS Optical Inserts are not suitable for driving or operating heavy machinery.

Do not use ZEISS Optical Inserts while driving, operating heavy machinery, or in any other situation that requires attention to safety. Disregard can cause serious accidents.

- Only use ZEISS Optical Inserts for their intended use.

2. ZEISS Optical Inserts must only be used when properly inserted into Apple Vision Pro.

Do not use ZEISS Optical Inserts for vision correction outside of Apple Vision Pro and / or as substitute for your prescription eyeglass lenses, contact lenses, or any other vision correction eyewear that is intended for normal wear and outside of Apple Vision Pro, as this could lead to unsafe situations. Disregard can cause serious accidents.


- Use ZEISS Optical Inserts only while wearing Apple Vision Pro.

3. ZEISS Optical Inserts are not intended for use as personal protective equipment.

ZEISS Optical Inserts are not classified as personal protective equipment. They are not intended for use as personal protective equipment. They do not claim to protect the eye from physical impact as they are neither shatterproof nor unbreakable. Disregarding intended use instructions may result in serious accidents.

- Only use ZEISS Optical Inserts for their intended use.

4. ZEISS Optical Inserts contain magnets that could impair certain medical devices and equipment.

 Improper use of ZEISS Optical Inserts could cause harm to the wearer or lead to unsafe situations.

- Always keep ZEISS Optical Inserts at least 20 mm away from an implanted electronic medical device (e.g., a pacemaker). Otherwise, the implanted electronic medical device could malfunction temporarily.

- If you have any questions, contact your health care professional for more information.

 **WARNING | Injury**

1. ZEISS Optical Inserts should be clean before use.


Proper use cannot be guaranteed if ZEISS Optical Inserts are not clean. This can lead to unsafe situations and product damage.

- Make sure ZEISS Optical Inserts are clean before use. Follow the maintenance and cleaning instructions in chapter 7 if cleaning is needed.

2. ZEISS Optical Inserts are not suitable to wear if broken or damaged.

If ZEISS Optical Inserts are broken or otherwise damaged, it could lead to unsafe situations (e.g., fragments could cause eye injuries or parts may be inadvertently swallowed).

- Do not wear ZEISS Optical Inserts in case they are broken or otherwise damaged.
- Keep ZEISS Optical Inserts out of the reach of children.
- Never swallow any parts of ZEISS Optical Inserts.

 **CAUTION | Potential risk or product damage**

1. ZEISS Optical Inserts *Prescription* are intended to be used only with your individual prescription.

Incorrect ZEISS Optical Inserts parameters (prescription and / or position of wear) can cause visual discomfort — eye strain, blurred vision, etc.

- If you experience any of these symptoms, stop using ZEISS Optical Inserts immediately. Contact your eye care professional to get your prescription checked.

2. ZEISS Optical Inserts are not to be modified or adjusted.

ZEISS Optical Inserts may not be modified as this could lead to product damage. In case of any modification of ZEISS Optical Inserts, the warranty expires.

- Do not modify ZEISS Optical Inserts.

5 Installation and use

Before the installation, read the ZEISS Optical Inserts instructions to ensure safe use and check the orientation marks on ZEISS Optical Inserts. With the help of the ZEISS Optical Inserts marked left (L) and right (R), ensure proper installation in Apple Vision Pro.

Get started:

1. Place each ZEISS Optical Insert into position over the displays. ZEISS Optical Inserts will snap in magnetically when properly attached over the Apple Vision Pro displays.
2. If prompted, scan the code included on the card to pair your Apple Vision Pro to the ZEISS Optical Inserts.
3. If necessary, the system may prompt you to swap to the © Light Seal Cushion which is provided separately in the Apple Vision Pro kit. The Light Seal Cushion is designed to maintain the appropriate fit of Apple Vision Pro. Consult the Apple Vision Pro User Guide for further instructions.

To safely remove the ZEISS Optical Inserts from your Apple Vision Pro, press your fingertip into the indentation on the side and gently lift the insert out.

7 Maintenance and cleaning



CAUTION | Product Damage

- 1. ZEISS Optical Inserts need to be fully dry before inserting into Apple Vision Pro.**

If ZEISS Optical Inserts are not clean, or are damp or wet, the proper use of ZEISS Optical Inserts cannot be guaranteed. This can lead to product damage.

 - Ensure ZEISS Optical Inserts are fully dry and clean before inserting into Apple Vision Pro.
- 2. Remove ZEISS Optical Inserts from Apple Vision Pro for cleaning.**

Do not clean ZEISS Optical Inserts while they are inserted in Apple Vision Pro. Improper cleaning can cause product damage.

 - Gently remove the ZEISS Optical Inserts from your Apple Vision Pro before cleaning.

4. Gently dry with dry side of microfiber cloth.
Ensure ZEISS Optical Inserts are fully dry before inserting into Apple Vision Pro.


Without ZEISS cleaning solutions

1. Rinse ZEISS Optical Inserts under lukewarm water.
 2. Use a pH-neutral cleaning agent.
 3. Gently dry with a clean microfiber cloth.
- Ensure ZEISS Optical Inserts are fully dry before inserting into Apple Vision Pro.

8 Troubleshooting

Issue(s):	<ul style="list-style-type: none"> • My vision is blurred, not sharp or unclear. • I perceive visual discomfort.
Try this:	<ul style="list-style-type: none"> • Double check if you are using ZEISS Optical Inserts with your correct eyeglass lens prescription. • Check if ZEISS Optical Inserts are clean, dry and free of dust and dirt. • Ensure that ZEISS Optical Inserts are installed properly. • Check if your eyeglass lens prescription may have changed. Please contact your eye care professional for more information.

11 Packaging label

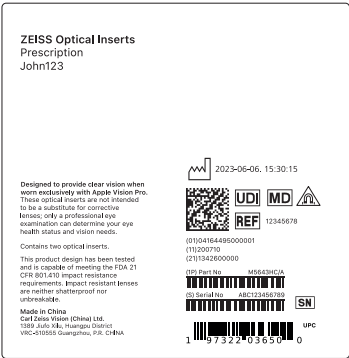

Carl Zeiss Vision International GmbH
Turnstrasse 27
73430 Aalen, Germany
zeiss.com/vision-care

6 Storage

When ZEISS Optical Inserts are not inserted in Apple Vision Pro, store ZEISS Optical Inserts safely to prevent scratches and other damage to the product.



CAUTION | Product Damage

- 1. ZEISS Optical Inserts are not suitable for extreme temperatures.**

Storing ZEISS Optical Inserts at temperatures above 45°C / 113°F (e.g., on the dashboard of a car in summer or in a sauna) can cause damage to ZEISS Optical Inserts.

 - Do not expose ZEISS Optical Inserts to extreme temperatures.
- 2. ZEISS Optical Inserts are not suitable for extreme temperature and humidity fluctuations.**

ZEISS Optical Inserts should not be exposed to extreme temperature and humidity fluctuations. This can cause damage to ZEISS Optical Inserts.

 - Do not expose ZEISS Optical Inserts to sudden temperature or humidity changes.

3. Use only described cleaning methods and agents to safely clean the ZEISS Optical Inserts.
- Improper cleaning can cause product damage.
- Do not use any cleaning methods and / or agents not described in this chapter.

ZEISS Lens Wipes

1. Lightly brush away dust and dirt with the folded lens wipe.
2. Unfold the lens wipe and gently clean the surface using circular movements until clean and dry.
3. Ensure ZEISS Optical Inserts are fully dry before inserting into Apple Vision Pro.

ZEISS Cleaning Spray

1. When heavily soiled, rinse ZEISS Optical Inserts under lukewarm water.
2. Apply ZEISS Cleaning Spray.
3. Use a clean microfiber cloth to gently clean the surface with circular movements.

9 Technical information

Refractive index	1.67
Coating	Functional multicoat
Optimization of power values	ZEISS Optical Inserts <i>Prescription</i> may be optimized from the original ordered prescription. In this case the refractive parameters may differ from the original order to provide the best visual experience. If you have any issues contact support.

10 Warranty

ZEISS Optical Inserts are provided according to applicable statutory and contractual stipulations in your country. Never modify ZEISS Optical Inserts. Any problems or damages to the ZEISS Optical Inserts caused by modifications or improper use will not be covered by the warranty.

12 Disposal and recycling

When this product reaches its end of life, dispose of it as instructed by your local authorities. According to local, state and federal law regulations, it is necessary to dispose of your product (packaging waste and ZEISS Optical Inserts) separately. Proper disposal will help conserve natural resources and protect both public health and the environment.

13 Contact information

Incident reporting and complaints:

For reporting an incident, causing injury from ZEISS Optical Inserts
contact: regulatoryaffairs.vision.us@zeiss.com

Legal Manufacturer

Carl Zeiss Vision International GmbH
Turnstrasse 27
73430 Aalen, Germany
zeiss.com/vision-care

For support and service view:

<https://support.apple.com/guide/apple-vision-pro>